

CareSens[®] PRO

Blood Glucose Monitoring System



*Actual Size

No Oxygen Interference

- Proven accuracy & precision
- Fresh capillary & venous whole blood



Specifications

- 0.4 μ L / 5 seconds
- 10-600 mg/dL (0.6-33.3 mmol/L)
- Hematocrit: 15-65%
- 1,000 test result memory
- 5-45°C (41-113°F)
- Fresh capillary / venous whole blood

Features

- Enzyme : GDH-FAD
- No coding required
- Strip ejection button
- Hypoglycemia indicator
- Strip expiration date indicator
- Auto back-lighted LCD (option)
- 4 alarms (PP2* and 3 time-set alarms)
*2 hours after meal
- 1, 7, 14, 30, 60, 90-day test averages
- Pre-meal, Post-meal, Fasting and No-mark flagging
- Data port

Easy Diabetes Care: Less work and less worry with CareSens® PRO

One-step Operation


- No-coding: Just insert the strip and you are ready to test in one step



① Insert test strip




② Apply blood sample

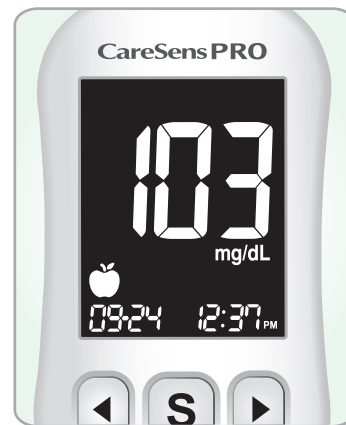

③ Result in 5 seconds

Improved Measurement Performance

- Advanced accuracy
 - 99.7% of results within ± 15 mg/dL (0.83 mmol/L) or 15% of the reference results, meeting ISO 15197:2013 (EN ISO 15197:2015) accuracy criteria
- Wide hematocrit range
 - 15-65%

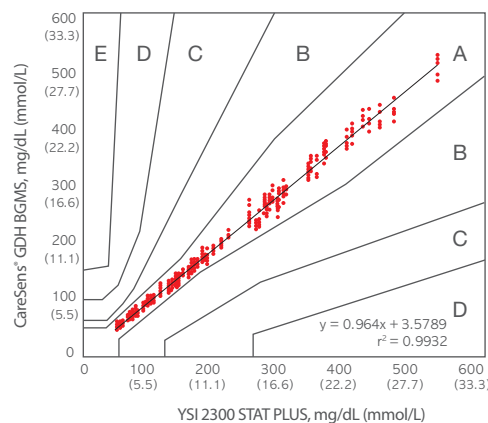

Practical Design

